

# Realism

## GLAZE PASTE

One For All-Ceramics



- Independent research and development by UPCERA
- Perfectly match with UPCERA zirconia and lithium disilicate blocks
- Easy cooperation with natural effect

### Glazing Program

Initial Temperature S °C	400-500
Preheating Time B mm:ss	4:00
Drying Time B mm:ss	2:00
Heating Rate T1 °C/mm	40-60
Holding Temperature T1 °C	760-830
Holding Time H1 mm	1:00
Cooling Rate T1 °C/mm	40-60

#### BASIC




#### INCISAL



#### FOSSA



#### EFFECT

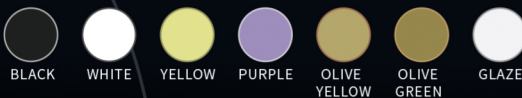

#### GINGIVAL



Shenzhen Upcera Dental Technology Co.,Ltd

Tel: +86-755-8882 0818

E-mail: sales@upcera.com

[www.upcera-dental.com](http://www.upcera-dental.com)

Follow us on Upcera

